




SKLARTECH 5000 / SKLARTECH 5000 MINIATURE

Laparoscopic Systems

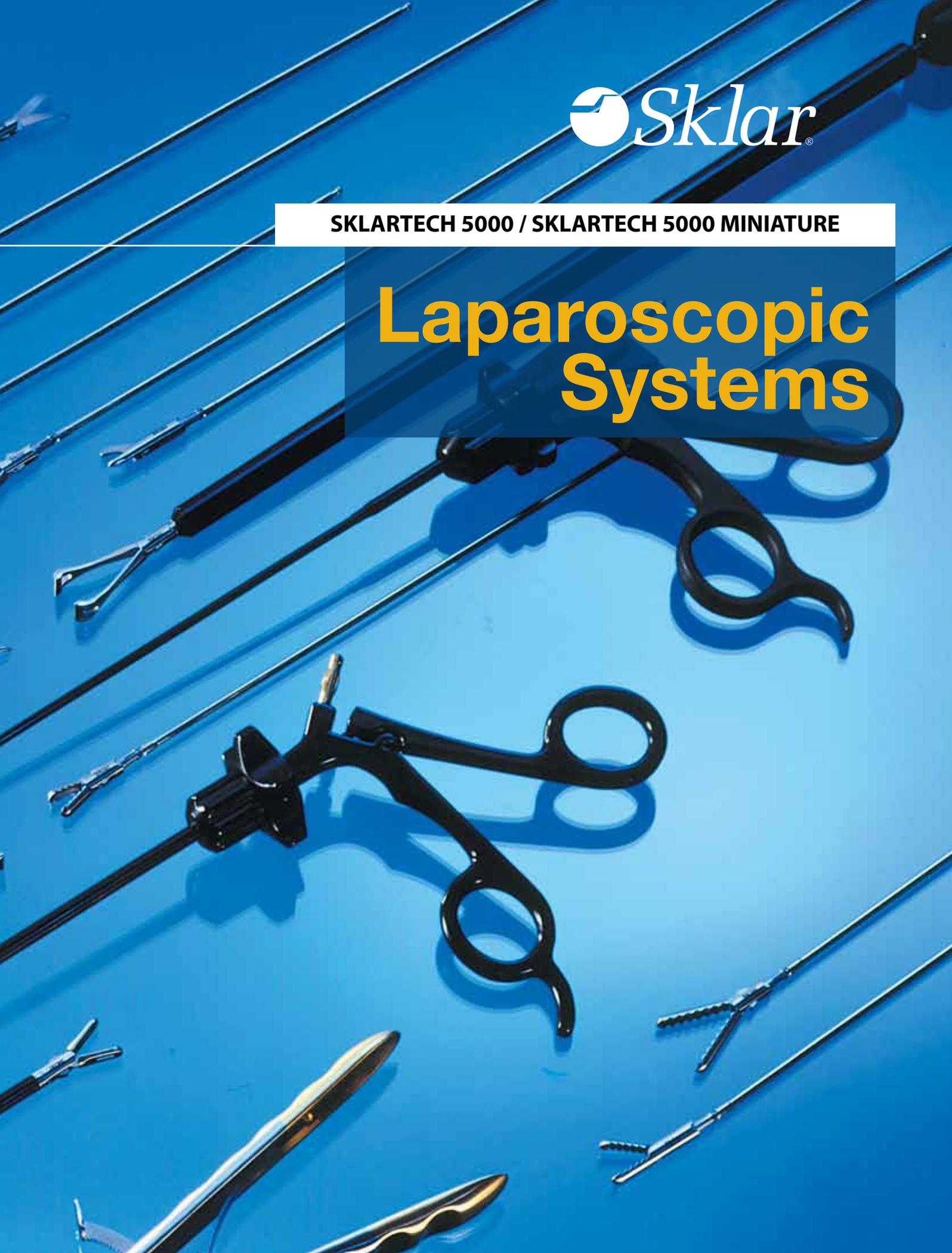

Sklartech 5000

Sklartech 5000

The Sklartech 5000 Laparoscopic Instrument line provides a full range of the most popular patterns. Instruments can be purchased as a complete set or individually. When purchased as a complete instrument, Sklartech Laparoscopic Scissors, Dissectors and Forceps include the insert with tip, insulated cannula and handle with locking collar.

As with all Sklar complete Laparoscopic instruments, the handles and replacement inserts are easily interchangeable to perform as intended and required by the surgeon. Each Sklartech 5000 scope is constructed of the highest quality German grade stainless steel to provide reliable durability and ease of use.

FEATURES:

INSULATED SHAFT provides vital protection against electrosurgical transfer.

INSERT WITH TIP is easily separated from cannula for cleaning and reprocessing with a convenient screw design.

LOCKING COLLAR features 360° rotation for ideal jaw position. Activate lock to provide tension and prevent accidental movement.

CARBON FIBER HANDLES are a comfortable ring construction that allows for precise hold of tissue. choose from ratchet (YC) or non-ratchet (XC) designs.

SINGLE PORT LAPAROSCOPY

Discover the **New Evolution** in minimally invasive Laparoscopic Surgery

Single Port Laparoscopic surgery is an innovative procedure in which the skilled surgeon operates almost exclusively through a single entry point, typically the patients navel. Unlike the traditional multi-port laparoscopic approach, SPL leaves only a single small scar and most patients report less discomfort and faster recovery.

- Decreased trauma and risk for herniation
- Enhanced cosmetic benefits with minimal visible scarring
- Reduced post-operative pain
- Accelerated recovery time



Item Number 31-9415
20mm Diameter

Trocar that supports (2) 5mm and (1) 10mm Trocar/Obturator

SINGLE PORT TROCARS

INSTRUMENT	DIAMETER	FEATURES
31-9415	20mm	Trocar that supports (2) 5mm and (1) 10mm Trocar/Obturator
31-9416	5mm	Trocar/Obturator Support
31-9417	10mm	Trocar/Obturator Support
31-9418	12.5mm	Trocar/Obturator Support




SUCTION IRRIGATION

INSTRUMENT	DIAMETER	LENGTH	FEATURES
31-9411	5mm	35cm	Non-Ratcheted Handle



ETH NEEDLE HOLDERS



INSTRUMENT	DIAMETER	LENGTH	FEATURES
31-9404	5mm	35cm	Single Curve, Non-Ratcheted Handle
31-9412	5mm	35cm	Double Curve, Non-Ratcheted Handle

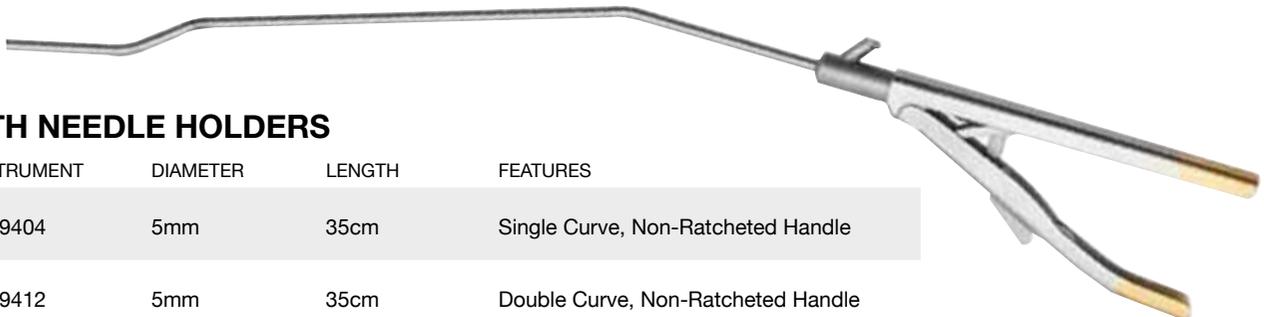



SINGLE CURVE INSTRUMENTS

INSTRUMENT	DIAMETER	LENGTH	FEATURES
31-9400	5mm	35cm	Metzenbaum Scissor, Non-Ratcheted Handle
31-9401	5mm	35cm	Maryland Dissector, Non-Ratcheted Handle
31-9402	5mm	35cm	Endo Clinch Grasping Forcep, Ratchet Handle
31-9403	5mm	35cm	Maxi Grasping Forcep, Ratchet Handle

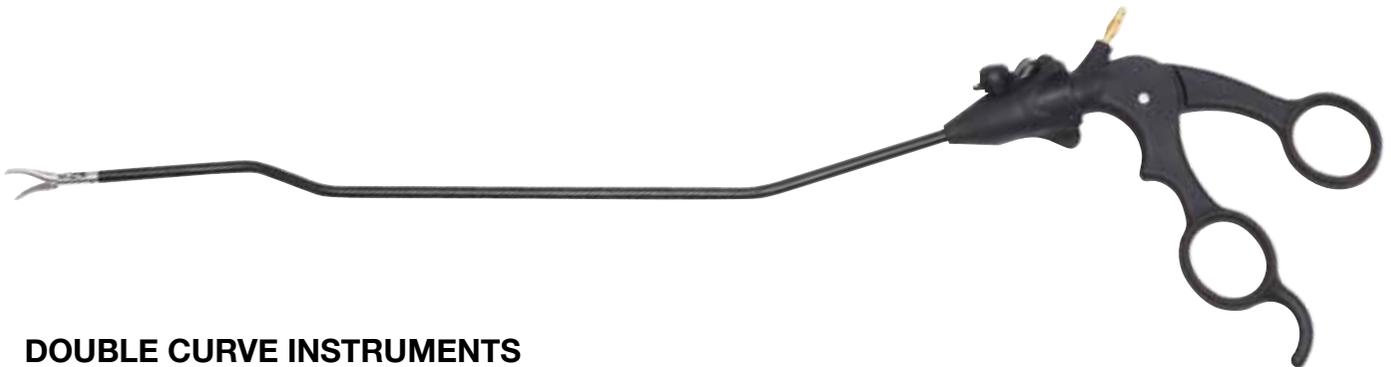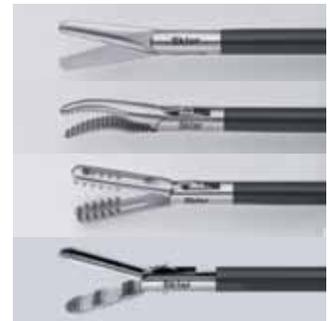

DOUBLE CURVE INSTRUMENTS

INSTRUMENT	DIAMETER	LENGTH	FEATURES
31-9405	5mm	35cm	Metzenbaum Scissor, Non-Ratcheted Handle
31-9406	5mm	35cm	Maryland Dissector, Non-Ratcheted Handle
31-9407	5mm	35cm	Endo Clinch Grasping Forcep, Ratcheted Handle
31-9408	5mm	35cm	Maxi Grasping Forcep, Ratcheted Handle
31-9409	5mm	35cm	L-Hook Electrode, Non-Ratcheted Handle
31-9410	5mm	35cm	Button Electrode, Non-Ratcheted Handle



Offering a full range of the most popular patterns and sizes, SklarTech 5000 Laparoscopic complete instruments come standard with insert and tip, insulated cannula and handle with locking collar. Sklar also carries a full line of interchangeable replacement parts, accessories and miniature instruments for your surgical preference.

PERITONEAL SCISSORS

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9026XC	31-9026	5mm	33cm	Straight, Blunt Tips, XC Handle



OPERATING SCISSORS

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9027XC	31-9027	10mm	33cm	Straight, XC Handle



METZENBAUM SCISSORS

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9005XC	31-9005	5mm	33cm	Straight, XC Handle
31-9005YC	31-9005	5mm	33cm	Straight, YC Handle



COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9010XC	31-9010	5mm	33cm	Curved, XC Handle
31-9010XC-45	31-9010-45	5mm	45cm	Curved, XC Handle



COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9012XC	31-9012	5mm	33cm	Curved, Serrated, XC Handle
31-9012XC-45	31-9012-45	5mm	45cm	Curved, Serrated, XC Handle



HOOK SCISSORS

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9016XC	31-9016	5mm	33cm	Standard, Flush Tips, XC Handle



COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9015XC	31-9015	5mm	33cm	Sharp Cutting Edge, Overlapping Tips, XC Handle
31-9015XC-45	31-9015-45	5mm	45cm	Sharp Cutting Edge, Overlapping Tips, XC Handle


MICRO SCISSORS

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9020XC	31-9020	5mm	33cm	Straight, Extremely Delicate, XC Handle



COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9025XC	31-9025	5mm	33cm	Curved, Extremely Delicate, XC Handle



MARYLAND DISSECTORS

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9030XC	31-9030	5mm	33cm	Cross Serration, XC Handle
31-9030XC-45	31-9030-45	5mm	45cm	Cross Serration, XC Handle



COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9031XC	31-9031	5mm	33cm	Diamond Serration, XC Handle
31-9120XC	31-9120	10mm	33cm	Diamond Serration, XC Handle


MIXTER FORCEPS

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9033XC	31-9033	5mm	33cm	90° Angle, Cross Serration, XC Handle



COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9035XC	31-9035	5mm	33cm	90° Angle, Diamond Serration XC Handle
31-9125XC	31-9125	10mm	33cm	90° Angle, Diamond Serration XC Handle
31-9125XC-45	31-9125-45	10mm	45cm	90° Angle, Diamond Serration XC Handle


DISSECTING FORCEP

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9036XC	31-9036	5mm	33cm	Curved, Serrated, XC Handle


DOLPHIN FORCEP

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9040XC	31-9040	5mm	33cm	Short, XC Handle



KELLY FORCEP

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9043XC	31-9043	5mm	33cm	Double Action, XC Handle


GLASSMANN FORCEP

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9096YC	31-9096	5mm	33cm	Atraumatic Clamp, YC Handle


TENACULUM FORCEPS

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9111YC	31-9111	5mm	33cm	1x1 Sharp Teeth, YC Handle
31-9113YC	31-9113	10mm	33cm	1x1 Sharp Teeth, YC Handle


KOCHER FORCEP

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9128XC	31-9128	10mm	33cm	60mm Jaws, Axial Serration, XC Handle


SPOON FORCEP

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9135XC	31-9135	10mm	33cm	XC Handle


CLAW FORCEPS

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9045XC	31-9045	5mm	33cm	Traumatic Jaws, XC Handle
31-9045YC	31-9045	5mm	33cm	Traumatic Jaws, YC Handle



COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9130XC	31-9130	10mm	33cm	Traumatic Jaws, XC Handle



BABCOCK FORCEPS

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9105YC	31-9105	5mm	33cm	YC Handle



COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9140YC	31-9140	10mm	33cm	YC Handle
31-9140XC	31-9140	10mm	33cm	XC Handle
31-9140YC-45	31-9140-45	10mm	45cm	YC Handle


ALLIS FORCEPS

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9067YC	31-9067	5mm	33cm	TC Inserted Jaws, Atraumatic Clamp, YC Handle



COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9144YC	31-9144	10mm	33cm	TC Inserted Jaws, Diamond Serrations, YC Handle


STRONG ATRAUMATIC CLAMP

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9145YC	31-9145	10mm	33cm	40mm Jaws, Straight, YC Handle


NEEDLE HOLDER

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9112YC	31-9112	5mm	33cm	10mm Jaws, Pointed, YC Handle



GRASPING FORCEPS

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9050XC	31-9050	5mm	33cm	Standard, Cross Serrations, XC Handle
31-9050YC	31-9050	5mm	33cm	Standard, Cross Serrations, YC Handle
31-9050YC-45	31-9050-45	5mm	45cm	Standard, Cross Serrations, YC Handle



COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9051YC	31-9051	5mm	33cm	Tapered, Rounded Jaws, YC Handle



COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9055YC	31-9055	5mm	33cm	Endo Clinch, Extremely Atraumatic Serrations, YC Handle



COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9080YC	31-9080	5mm	33cm	With Spoon, YC Handle
31-9080XC	31-9080	5mm	33cm	With Spoon, XC Handle



COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9085YC	31-9085	5mm	33cm	Pointed with Spoon, YC Handle



COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9100YC	31-9100	5mm	33cm	Maxi (Fundus), YC Handle



WAVE FORCEPS

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9060YC	31-9060	5mm	33cm	Atraumatic Wave like Serrations, Single Action, YC Handle
31-9060YC-45	31-9060-45	5mm	45cm	Atraumatic Wave like Serrations, Single Action, YC Handle
31-9061YC	31-9061	5mm	33cm	Atraumatic Wave like Serrations, Double Action, YC Handle


DEBAKEY FORCEPS

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9065YC	31-9065	5mm	33cm	40mm Jaws, YC Handle
31-9065XC	31-9065	5mm	33cm	40mm Jaws, XC Handle
31-9065YC-45	31-9065-45	5mm	45cm	40mm Jaws, YC Handle
31-9065XC-45	31-9065-45	5mm	45cm	40mm Jaws, XC Handle


BABCOCK/DEBAKEY FORCEP

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9142XC	31-9142	10mm	33cm	Atraumatic Serrations, XC Handle


COBRA FORCEP

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9075YC	31-9075	5mm	33cm	2x3 Teeth, Double Action, YC Handle


ALLIGATOR FORCEP

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9070YC	31-9070	5mm	33cm	Double Action, YC Handle



ATRAUMATIC GRASPERS

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9098YC	31-9098	5mm	33cm	Single Action, YC Handle



COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9097XC	31-9097	5mm	33cm	Diamond Serrations, XC Handle
31-9097YC	31-9097	5mm	33cm	Diamond Serrations, YC Handle


FENESTRATED GRASPING DISSECTOR

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9042XC	31-9042	5mm	33cm	16mm Double Action Jaw, XC Handle


FENESTRATED GRASPING FORCEPS

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9090YC	31-9090	5mm	33cm	22mm Double Action Jaw, YC Handle
31-9090XC-45	31-9090-45	5mm	45cm	22mm Double Action Jaw, XC Handle



COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9095XC	31-9095	5mm	33cm	40mm Double Action Jaw, XC Handle
31-9095YC	31-9095	5mm	33cm	40mm Double Action Jaw, YC Handle
31-9095XC-45	31-9095-45	5mm	45cm	40mm Double Action Jaw, XC Handle
31-9095YC-45	31-9095-45	5mm	45cm	40mm Double Action Jaw, YC Handle



PUNCH

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9046XC	31-9046	5mm	33cm	Probe Excision


BIOPSY PUNCH

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9047XC	31-9047	5mm	33cm	No Tooth, XC Handle



COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9110YC	31-9110	5mm	33cm	With Tooth, YC Handle
31-9110YC-45	31-9110-45	5mm	45cm	With Tooth, YC Handle




COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-9115XC	31-9115	5mm	33cm	Spoon with Teeth, XC Handle
31-9115YC	31-9115	5mm	33cm	Spoon with Teeth, YC Handle
31-9115YC-45	31-9115-45	5mm	45cm	Spoon with Teeth, YC Handle



SKLARLITE™ Endoscopy Rack

With Silicone Holding Clamps

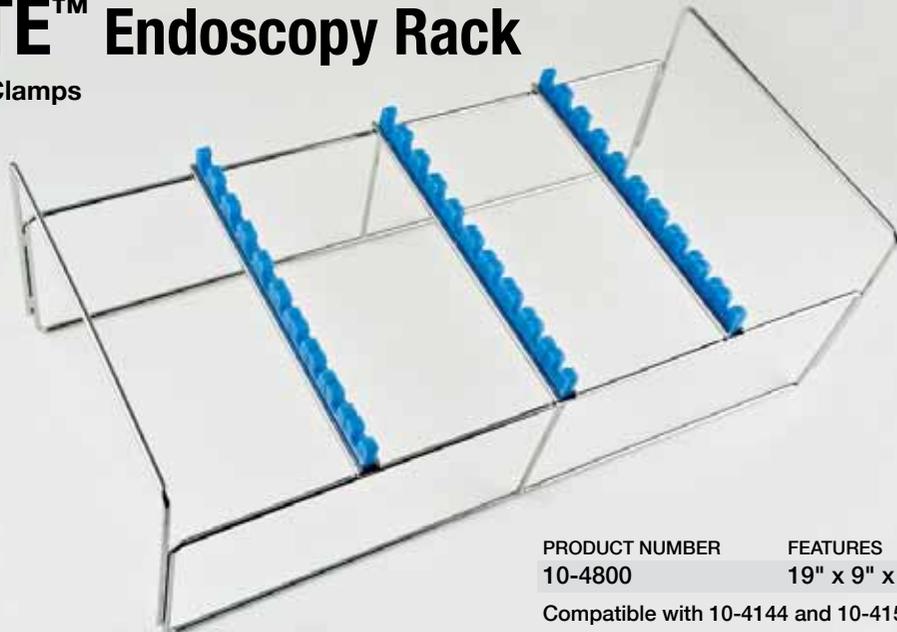

PRODUCT NUMBER	FEATURES
10-4800	19" x 9" x 6½"
Compatible with 10-4144 and 10-4156	

SKLAR “XC” CARBON HANDLE

INSTRUMENT	FEATURES
31-9001XC	Standard, non-ratchet, lightweight, carbon fiber reinforced non-conductive safety handle with locking collar. Fully rotating with a proximally positioned high-frequency connector.

SKLAR “YC” CARBON HANDLE

INSTRUMENT	FEATURES
31-9002YC	Standard, ratchet, lightweight, carbon fiber reinforced non-conductive safety handle with locking collar. Fully rotating with a proximally positioned high-frequency connector.

REPLACEMENT LOCKING COLLAR

INSTRUMENT	FEATURES
31-9160	For Sklar Standard Handles



Instruments may be ordered as complete sets or assembled from components ordered separately. Optional handles for customizing instruments – interchangeable for 3mm, 5mm or 10mm – are available for special order. Contact your Sklar Account Manager for details.

CANNULAS

INSTRUMENT	DIAMETER	LENGTH	FEATURES
31-9155	5mm	33cm	Insulated
31-9155-45	5mm	45cm	Insulated



31-9150	10mm	33cm	Insulated With Permanent Locking Collar Attached
31-9150-45	10mm	45cm	Insulated With Permanent Locking Collar Attached



HASSON TROCAR WITH TRUMPET VALVE

COMPLETE INSTRUMENT	REPLACEMENT OBTURATOR	LENGTH	FEATURES
31-1002	31-0004	11mm	Hasson Trumpet Valve Cannula with 11mm Blunt Tip Obturator


STANDARD TROCAR WITH TRUMPET VALVE

COMPLETE INSTRUMENT	REPLACEMENT OBTURATOR	LENGTH	FEATURES
31-1007	31-1010	5.5mm	Trumpet Valve Cannula With 5.5mm Sharp Obturator


STANDARD TROCAR WITH TRAP DOOR

COMPLETE INSTRUMENT	REPLACEMENT OBTURATOR	LENGTH	FEATURES
31-1004	31-1011	11mm	Trap Door Cannula With 11mm Sharp Obturator


STANDARD TROCAR WITH AUTOMATIC VALVE

COMPLETE INSTRUMENT	REPLACEMENT OBTURATOR	LENGTH	FEATURES
31-1015	31-1010	5.5mm	Automatic Valve Cannula with 5.5mm Sharp Obturator, Dismountable
31-1016	31-1011	11mm	Automatic Valve Cannula with 11mm Sharp Obturator, Dismountable



STANDARD TROCAR WITH SILICON VALVE

COMPLETE INSTRUMENT	REPLACEMENT OBTURATOR	LENGTH	FEATURES
31-1017	31-1010	5.5mm	Silicon Valve Cannula with LL Stopcock & 5.5mm Sharp Obturator, Dismountable
31-1018	31-1011	11mm	Silicon Valve Cannula with LL Stopcock & 11mm Sharp Obturator, Dismountable



OBTURATORS

INSTRUMENT	LENGTH	FEATURES
31-0004	11mm	Blunt Tip Obturator for 11mm Hasson Cannula
31-1011	11mm	Extremely Sharp Pyramidal Tip Obturator for Standard Trocar
31-1010	5.5mm	Extremely Sharp Pyramidal Tip Obturator for Standard Trocar

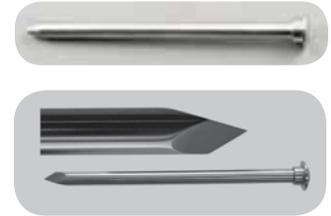
THREADED CANNULA SLEEVE

INSTRUMENT	LENGTH	FEATURES
31-1020	5.5mm	Interchangeable with Standard Smooth Cannulas


SNAP ON FLIP REDUCER

INSTRUMENT	LENGTH	FEATURES
31-1042	-	Reduces from 11mm to 5.5mm


S-RETRACTOR

INSTRUMENT	LENGTH	FEATURES
60-6130	6¼"	Set of 2


TROCAR CANNULA PARTS

FOR USE WITH STANDARD & HASSON VALVE CANNULAS				
EXTERNAL SEALS	5.5mm	11mm	11mm (Hasson)	11mm (Hasson)
	RED	BLUE	BLUE	GRAY
	31-0003	31-0001	31-0001	31-0002
FOR USE WITH TRAP DOOR CANNULAS				
INTERNAL SEAL	5.5mm		11mm	
	GRAY		YELLOW	
	31-4323		31-4324	
FOR USE WITH SILICON VALVE CANNULAS				
SILICON VALVE	5.5mm		11mm	
	31-4328		31-4329	




**SKLARTECH
5000**

Each coagulating Electrode is able to connect to the suction control valve of your choice (Listed below). Insulated Electrodes and Control Valves sold separately.


INSULATED ELECTRODES, WITH SUCTION

INSTRUMENT	DIAMETER	FEATURES
31-2235 (A)	5mm	Spatula
31-2237 (B)	5mm	L-Hook
31-2241 (C)	5mm	Knife
31-2243 (D)	5mm	J-Hook
31-2244 (E)	5mm	Blunt

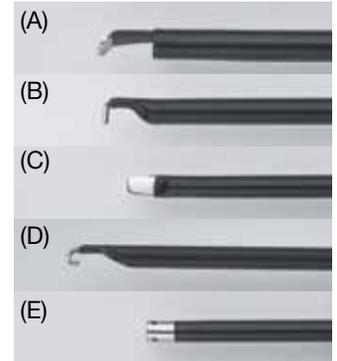

HIGH FREQUENCY SUCTION CONTROL VALVES


INSTRUMENT	DIAMETER	FEATURES
31-2227 (A)	5mm	Insulated Suction Control with Single Trumpet Valve
31-2233 (B)	5mm	Insulated Suction Control with Double Trumpet Valve
31-2229 (C)	5mm	Insulated Suction Control with Double Trumpet Valve and Handle

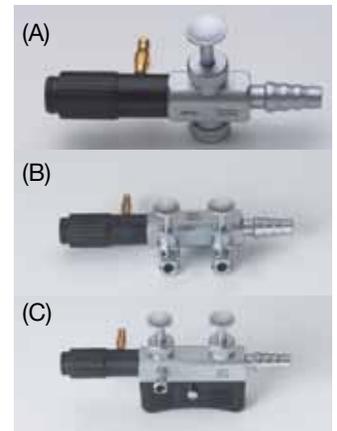

INSULATED ELECTRODES, WITHOUT SUCTION

Each Electrode has a permanently fitted adapter for connecting to the Electrosurgical Cable. Insulated Electrodes and Cables sold separately.

INSTRUMENT	DIAMETER	FEATURES
31-2435 (A)	5mm	Spatula
31-2437 (B)	5mm	L-Hook
31-2443 (C)	5mm	J-Hook
31-2445 (D)	5mm	Button
31-2447 (E)	5mm	Cone

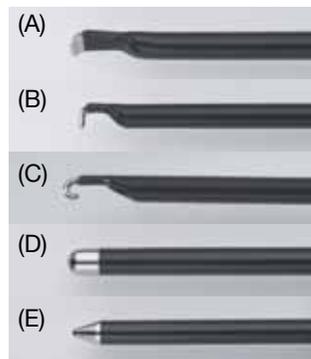

ELECTROSURGICAL CABLE

INSTRUMENT	FEATURES
31-2242	9' 4mm Tip, Monopolar: Compatible with Martin, Berchtold, Erbe T-Series
31-2253	10' 8mm Tip, Monopolar: Compatible with ValleyLab, Bovie, Conmed



METZENBAUM SCISSORS

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-4306XC	31-4306	3mm	33cm	Mini, Straight, XC Handle
31-4303XC	31-4303	3mm	33cm	Mini, Curved, XC Handle
31-4307XC	31-4307	3mm	33cm	Mini, Curved, Serrated, XC Handle


HOOK SCISSORS

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-4304XC	31-4304	3mm	33cm	Mini, Delicate, Straight, XC Handle


MARYLAND DISSECTOR

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-4310XC	31-4310	3mm	33cm	Mini, Cross Serrations, Curved XC Handle



FORCEPS

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-4314YC	31-4314	3mm	33cm	Mini Duck Nose, Cross Serrations, Straight, YC Handle
31-4309XC	31-4309	3mm	33cm	Mini Kelly, Double Action Delicate Jaws, XC Handle
31-4315XC	31-4315	3mm	33cm	Mini Spoon, XC Handle
31-4316XC	31-4316	3mm	33cm	Mini Spoon, Serrated, XC Handle



GRASPING FORCEPS

COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-4308YC	31-4308	3mm	33cm	Mini Grasping Forcep, YC Handle
31-4311XC	31-4311	3mm	33cm	Mini Atraumatic, Dissector, Straight XC Handle
31-4312YC	31-4312	3mm	33cm	Mini Grasping Forcep, Heavy Serrations YC Handle
31-4313YC	31-4313	3mm	33cm	Mini Cobra, Heavy Serrations YC Handle
31-4317YC	31-4317	3mm	33cm	Mini Fenestrated, 8mm Double Action Jaws, YC Handle

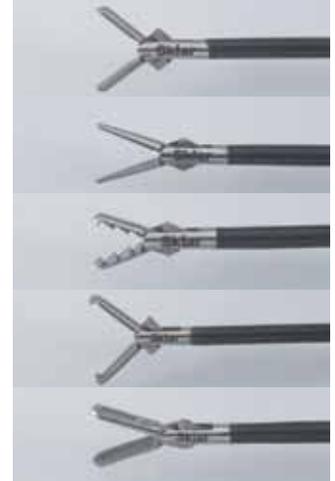

PUNCH


COMPLETE INSTRUMENT	REPLACEMENT INSERT	DIAMETER	LENGTH	FEATURES
31-4305XC	31-4305	3mm	33cm	Mini Punch, Delicate



NEEDLE HOLDER

NON-INSULATED WITH TUNGSTEN CARBIDE JAW INSERTS AND AXIAL HANDLE

COMPLETE INSTRUMENT	DIAMETER	LENGTH	FEATURES
31-9360AX	3mm	33cm	Mini, Straight, Axial Ring Handle
31-9362AX	3mm	33cm	Mini, Curved Right, Axial Ring Handle
31-9364AX	3mm	33cm	Mini, Curved Left, Axial Ring Handle

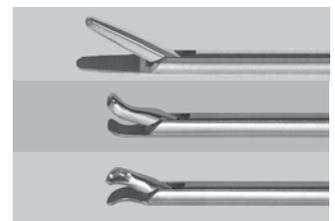

MINIATURE ENDOSCOPE

COMPLETE INSTRUMENT	DIAMETER	LENGTH	FEATURES
31-4301	2.7mm	30cm	Miniature, Autoclavable



MINIATURE TROCAR

COMPLETE INSTRUMENT	REPLACEMENT OBTURATOR	DIAMETER	FEATURES
31-4302	31-4318	3mm	Mini, with Trap Door Cannula, Reusable, Includes 3mm Obturator
31-4320	-	3mm	External Seal, Red For use with Standard and Hasson Cannulas
31-4321	-	3mm	Internal Seal, Gray For use with Trap Door Cannulas only




ELECTROSURGERY

INSULATED ELECTRODES WITHOUT SUCTION

Each Electrode has a permanently fitted adapter for connecting to the Electrosurgical Cable. Insulated Electrodes and Cables sold separately.

INSTRUMENT	DIAMETER	FEATURES
31-4325 (A)	3mm	Mini, Spatula
31-4326 (B)	3mm	Mini, L-Hook
31-4327 (C)	3mm	Mini, J-Hook

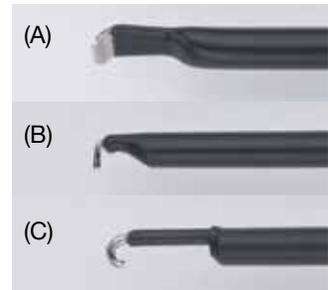

ELECTROSURGICAL CABLE

INSTRUMENT	LENGTH	FEATURES
31-2242	9'	4mm Tip, Monopolar: Compatible with Martin, Berchtold, Erb T-Series
31-2253	10'	8mm Tip, Monopolar: Compatible with ValleyLab, Bovie, Conmed





VERRES NEEDLE WITH STOPCOCK

INSTRUMENT	DIAMETER	LENGTH	FEATURES
31-2327	2.3mm	100mm	Stopcock
31-2350	2.3mm	120mm	Stopcock
31-2352	2.3mm	150mm	Stopcock



NEEDLE FOR INJECTION AND IRRIGATION

INSTRUMENT	DIAMETER	LENGTH	FEATURES
31-2355	5mm	33cm	Luer Lock, 2mm Hollow Needle
31-2355-45	5mm	45cm	Luer Lock, 2mm Hollow Needle



GRADUATED PROBE

INSTRUMENT	DIAMETER	LENGTH	FEATURES
31-2365	5.5mm	33cm	Graduated



ESOPHAGUS RETRACTOR

INSTRUMENT	DIAMETER	LENGTH	FEATURES
31-4565	5mm	33cm	For Gastric Band Application



FAN RETRACTOR

INSTRUMENT	DIAMETER	LENGTH	FEATURES
31-4562	5mm	33cm	3 Finger Articulating Fan, Slightly Bent



FAN RETRACTOR

INSTRUMENT	DIAMETER	LENGTH	FEATURES
31-4563	10mm	33cm	5 Finger Articulating Fan, Slightly Bent



SUCTION IRRIGATION SYSTEM

COMPLETE INSTRUMENT	DIAMETER	LENGTH	FEATURES
31-2238	5mm	33cm	Dual Trumpet Valve (also fits 10mm Sheath)
31-2238-45	5mm	45cm	Dual Trumpet Valve (also fits 10mm Sheath)
31-2258	5mm	33cm	Single Trumpet Valve (also fits 10mm Sheath)
31-4003	5mm	33cm	Sheath Only
31-4004	10mm	33cm	Sheath Only
31-2205	5mm	33cm	Tube Channel, for Laser-Fiber
31-2238OR	-	-	O-Ring, Replacement Pk/10



HIGH-VOLUME SUCTION IRRIGATION SYSTEM

COMPLETE INSTRUMENT	DIAMETER	LENGTH	FEATURES
31-2239	5mm 10mm	33cm	Dual Trumpet Valve (includes both 5mm and 10mm Sheaths)
31-2347	5mm	33cm	Sheath Only
31-2349	10mm	33cm	Sheath Only



TC NEEDLEHOLDER

NON-INSULATED WITH TUNGSTEN CARBIDE JAW INSERTS AND AXIAL HANDLE

COMPLETE INSTRUMENT	DIAMETER	LENGTH	FEATURES
31-2550	5mm	33cm	Standard, Straight, Flushable Shaft, 10mm Jaw Axial Handle
31-2550-45	5mm	45cm	Standard, Straight, Flushable Shaft, 10mm Jaw Axial Handle
31-2551	5mm	33cm	Standard, Curved Right, Flushable Shaft, 10mm Jaw Axial Handle
31-2552	5mm	33cm	Standard, Curved Left, Flushable Shaft, 10mm Jaw Axial Handle
31-2552-45	5mm	45cm	Standard, Curved Left, Flushable Shaft, 10mm Jaw Axial Handle
31-2553	5mm	33cm	Standard, Straight, Solid Shaft, 8mm Jaw Axial Handle



NEW

LIGATION CLIPS

Sklar is pleased to offer the most atraumatic shaped clip on the market. Sklar Ligation Clips' unique "U-Shape" design provides a stronger, more permanent hold. The unique and highly-advanced cross-serrations on each clip make it extremely secure, eliminating slippage from the vessel after closure. And the solid titanium gives a scatter-free CT image and is MRI compatible.

Sklar Ligation Clips are conveniently packaged in sterile cartridges holding six single clips, eliminating waste. Cartridges feature an adhesive backing to easily affix to any surface.

Perfect compatibility between special applicators and clips allow for an effortless loading system while the chevron shape facilitates a surgeon's choice "tip-first" closing action. Vessels are contained within the clip so full control is available before ligation.



STERILE

INSTRUMENT	SIZE	FEATURES	QUANTITIES
31-7300	Small	Yellow	6/pk, 20pk/cs, 120/cs
31-7301	Medium	Blue	6/pk, 20pk/cs, 120/cs
31-7302	Medium/Large	Green	6/pk, 20pk/cs, 120/cs
31-7303	Large	Orange	6/pk, 20pk/cs, 120/cs



OPEN SURGERY APPLIERS

- Easy and Convenient Loading System
- Spring Action Controls Jaw Compression
- Firm Grasp of Clips
- Safe and Precise Closing Action
- Color Coded for Easy Identification
- Full Range of Sizes to Fit any Application

SMALL

INSTRUMENT	LENGTH	FEATURES
31-7305	6" 150mm	Straight, Yellow
31-7306	7½" 190mm	Straight, Yellow
31-7307	7½" 190mm	Right Angle, Yellow

MEDIUM

INSTRUMENT	LENGTH	FEATURES
31-7308	6" 150mm	Straight, Blue
31-7309	7½" 190mm	Straight, Blue
31-7310	7½" 190mm	Right Angle, Blue
31-7311	11" 280mm	Straight, Blue
31-7312	11" 280mm	Right Angle, Blue

MEDIUM/LARGE

INSTRUMENT	LENGTH	FEATURES
31-7313	7½" 190mm	Straight, Green
31-7314	7½" 190mm	Right Angle, Green
31-7315	11" 280mm	Straight, Green
31-7316	11" 280mm	Right Angle, Green

LARGE

INSTRUMENT	LENGTH	FEATURES
31-7317	7½" 190mm	Straight, Orange
31-7318	7½" 190mm	Right Angle, Orange
31-7319	11" 280mm	Straight, Orange
31-7320	11" 280mm	Right Angle, Orange



Clips and Appliers sold separately.

ROTATABLE APPLIERS

- Easy and Convenient Loading System
- Full 360° Rotation of applicator shape
- Sized to match popular clip sizes


INSTRUMENT	DIAMETER	LENGTH	FEATURES
31-7321	5mm	35cm	Clip Applicator for Small Ligation Clips



INSTRUMENT	DIAMETER	LENGTH	FEATURES
31-7322	8mm	35cm	Clip Applicator for Medium Ligation Clips



INSTRUMENT	DIAMETER	LENGTH	FEATURES
31-7323	10mm	35cm	Clip Applicator for Medium/Large Ligation Clips

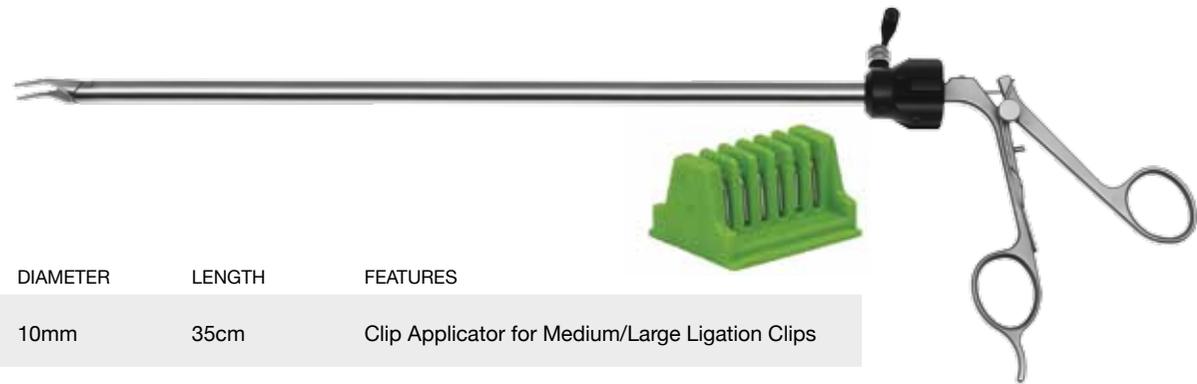

INSTRUMENT	DIAMETER	LENGTH	FEATURES
31-7324	12.5mm	35cm	Clip Applicator for Large Ligation Clips



Clips and Appliers sold separately.

HX-LINE Disposable Sheath/Insert with Reusable Handles

Sklar is proud to introduce the SKLARTECH 5000HX-Line of Disposable Laparoscopic Instruments featuring a disposable sheath/insert combination and reusable carbon fiber handle. With a disposable sheath/insert combination, all parts of the instrument that had direct contact during surgery can be disposed of immediately after procedures. The carbon fiber handles, available with or without ratchet, can quickly and easily be cleaned and sterilized for continued use.

- Choose from a full line of the most popular patterns.
- Benefit from factory-honed instrument edge with every use.
- Improved operating room efficiency.
- Eliminate difficult reprocessing procedures.
- Enhanced efficiency and cost savings.

NEW!

Sklartech 5000HX



REUSABLE HANDLES

INSTRUMENT	DESCRIPTION	QUANTITIES
31-7001XC	Carbon Fiber without Ratchet	Each
31-7002YC	Carbon Fiber with Ratchet	Each

DISPOSABLE SHEATH/INSERT

INSTRUMENT	DIAMETER	LENGTH	FEATURES	QUANTITIES
31-7070	5mm	33cm	Fenestrated Grasping Forcep, Serrated	10/pk
31-7070-45	5mm	45cm	Fenestrated Grasping Forcep, Serrated	10/pk
31-7030	5mm	33cm	Maryland Dissector, Curved, Cross Serrations	10/pk
31-7030-45	5mm	45cm	Maryland Dissector, Curved, Cross Serrations	10/pk
31-7055	5mm	33cm	Endo Clinch Forcep, Extremely Atraumatic, Serrated	10/pk
31-7055-45	5mm	45cm	Endo Clinch Forcep, Extremely Atraumatic, Serrated	10/pk
31-7061	5mm	33cm	Wave Grasper Forcep, Wavelike Serrations	10/pk
31-7061-45	5mm	45cm	Wave Grasper Forcep, Wavelike Serrations	10/pk
31-7035	5mm	33cm	Mixer Clamp Forcep, 90° Angle	10/pk
31-7035-45	5mm	45cm	Mixer Clamp Forcep, 90° Angle	10/pk
31-7040	5mm	33cm	Dolphin Nose Forcep, Serrated	10/pk
31-7040-45	5mm	45cm	Dolphin Nose Forcep, Serrated	10/pk
31-7073	5mm	33cm	Allis Forcep, 3x4 Teeth	10/pk
31-7073-45	5mm	45cm	Allis Forcep, 3x4 teeth	10/pk
31-7016	5mm	33cm	Hook Scissors, Single Action	10/pk
31-7016-45	5mm	45cm	Hook Scissors, Single Action	10/pk
31-7012	5mm	33cm	Metzenbaum Scissors, Curved, Serrated	10/pk
31-7012-45	5mm	45cm	Metzenbaum Scissors, Curved, Serrated	10/pk
31-7007	5mm	33cm	Metzenbaum Scissors, Straight, Serrated	10/pk
31-7007-45	5mm	45cm	Metzenbaum Scissors, Straight, Serrated	10/pk
31-7107	5mm	33cm	Babcock-Debaquey Forcep, Atraumatic Serrations	10/pk
31-7107-45	5mm	45cm	Babcock-Debaquey Forcep, Atraumatic Serrations	10/pk



Item Number 31-7030
Maryland Forcep, Curved, Cross Serrations

SKLARTECH 5000X

Sklar is pleased to introduce a new line of sterile, 100% Disposable Laparoscopic instruments. Complete with HF Port, the X-Line delivers quality and performance at an economical price. Improve safety and reduce the risk of cross contamination to patient and staff by disposing of the entire instrument immediately following procedures.

- Advantage of factory-hone instruments with every use.
- Simplify cost-accounting per operation.
- Eliminate reprocessing time and expense.
- Eliminate Maintenance.
- Reduce inventory.

NEW!

Sklartech 5000X



X-LINE DISPOSABLE INSTRUMENTS, COMPLETE UNIT

COMPLETE INSTRUMENT	DIAMETER	LENGTH	FEATURES	QUANTITIES
31-7012XC	5mm	33cm	Metzenbaum Scissors, Curved with Non-Ratcheted Handle, Red	6/pk
31-7030XC	5mm	33cm	Maryland Dissector, Cross Serrations with Non-Ratcheted Handle, Yellow	6/pk
31-7055XC	5mm	33cm	Endo Clinch Forcep, Cross Serrations with Non-Ratcheted Handle, Blue	6/pk
31-7055YC	5mm	33cm	Endo Clinch Forcep, Cross Serrations with Ratchet Handle, Blue	6/pk

DISPOSABLE TROCARS

INSTRUMENT	DIAMETER	LENGTH	FEATURES	QUANTITIES
31-7200	10mm	10cm	Hasson, Blunt	10/pk
31-7205	12.5mm	10cm	Hasson, Blunt	10/pk
31-7210	5mm	10cm	Safety Shield	10/pk
31-7215	10mm	10cm	Safety Shield	10/pk
31-7220	12.5mm	10cm	Safety Shield	10/pk
31-7225	5mm	10cm	Conical Tip, Pyramidal Obturator	10/pk
31-7230	10mm	10cm	Conical Tip, Conical Obturator	10/pk
31-7235	5mm	10cm	Dilating Tip Obturator	10/pk
31-7240	10mm	10cm	Dilating Tip Obturator	10/pk
31-7245	12.5mm	10cm	Dilating Tip Obturator	10/pk
31-7250	5mm	10cm	Sheath Only	10/pk
31-7255	10mm	10cm	Sheath Only	10/pk
31-7260	12.5mm	10cm	Sheath Only	10/pk



SKLARLITE™

STERILIZATION CONTAINER SYSTEM

Pre-Vacuum Steam Sterilization

The SklarLite™ Rigid Sterilization Container System provides an organized, routine and secure approach to the sterilization of surgical instruments. Customize your SklarLite™ System by selecting colored lids and ID labels. Available in your choice of six vibrant colors for easy identification of instrument sets by procedure and/or department. Choose from various configurations to fit your surgery center needs.

- Anodized aluminum alloy is strong and durable yet lightweight and easy to handle.
- Stackable Containers maximizes space in the autoclave and storage shelf space; reduces cycle runs.
- Stackable Baskets provide ease of instrument transfer per procedure and maximizes space inside containers.
- Silicon mats secure delicate and laparoscopic instruments to minimize damage.
- Tamper evident seals prevent accidental sterile field breach.
- Long Lasting Sterilization.



STANDARD CONTAINERS



Consists of a body and perforated lid and used with Sklar's disposable filters and other compatible accessories for sterilization. Choose from either perforated or non-perforated bottoms. Lids available in six different colors for easy identification by procedure or department.

SAFE CONTAINERS



Consists of identical body and perforated lid, but features an additional solid cover that safeguards the perforation field and contents during transportation and storage. Choose from either perforated or non-perforated bottoms. Safety lids available in six different colors for easy identification by procedure or department.

SKLARLITE™ HALF-SIZE CONTAINERS

11¼" L X 11" W

STANDARD CONTAINERS			SAFE CONTAINERS	
Non-Perforated Bottom	Perforated Bottom	Size	Non-Perforated Bottom	Perforated Bottom
10-4500	10-4512	4"	10-4524	10-4536
10-4502	10-4514	4"	10-4526	10-4538
10-4504	10-4516	4"	10-4528	10-4540
10-4506	10-4518	4"	10-4530	10-4542
10-4508	10-4520	4"	10-4532	10-4544
10-4510	10-4522	4"	10-4534	10-4546
10-4548	10-4560	5"	10-4572	10-4584
10-4550	10-4562	5"	10-4574	10-4586
10-4552	10-4564	5"	10-4576	10-4588
10-4554	10-4566	5"	10-4578	10-4590
10-4556	10-4568	5"	10-4580	10-4592
10-4558	10-4570	5"	10-4582	10-4594
10-4596	10-4608	6"	10-4620	10-4632
10-4598	10-4610	6"	10-4622	10-4634
10-4600	10-4612	6"	10-4624	10-4636
10-4602	10-4614	6"	10-4626	10-4638
10-4604	10-4616	6"	10-4628	10-4640
10-4606	10-4618	6"	10-4630	10-4642
10-4644	10-4656	8"	10-4668	10-4680
10-4646	10-4658	8"	10-4670	10-4682
10-4648	10-4660	8"	10-4672	10-4684
10-4650	10-4662	8"	10-4674	10-4686
10-4652	10-4664	8"	10-4676	10-4688
10-4654	10-4666	8"	10-4678	10-4690
10-4692	10-4704	10"	10-4716	10-4728
10-4694	10-4706	10"	10-4718	10-4730
10-4696	10-4708	10"	10-4720	10-4732
10-4698	10-4710	10"	10-4722	10-4734
10-4700	10-4712	10"	10-4724	10-4736
10-4702	10-4714	10"	10-4726	10-4738



ACCESSORIES

PRODUCT NUMBER	MODEL	FEATURES
10-4784	Half Size Wire Basket	9½" x 10" x 1"
10-4786	Half Size Wire Basket	9½" x 10" x 2"
10-4788	Half Size Wire Basket	9½" x 10" x 3"
10-4790	Half Size Wire Basket	9½" x 10" x 4"
10-4970	Half Size Silicon Mat	8½" x 9"



SKLARLITE™ MID-SIZE CONTAINERS
18¼" L X 11" W

STANDARD CONTAINERS			SAFE CONTAINERS	
Non-Perforated Bottom	Perforated Bottom	Size	Non-Perforated Bottom	Perforated Bottom
10-4300	10-4312	4"	10-4324	10-4336
10-4302	10-4314	4"	10-4326	10-4338
10-4304	10-4316	4"	10-4328	10-4340
10-4306	10-4318	4"	10-4330	10-4342
10-4308	10-4320	4"	10-4332	10-4344
10-4310	10-4322	4"	10-4334	10-4346
10-4348	10-4360	5"	10-4370	10-4382
10-4350	10-4361	5"	10-4372	10-4384
10-4352	10-4362	5"	10-4374	10-4386
10-4354	10-4364	5"	10-4376	10-4388
10-4356	10-4366	5"	10-4378	10-4390
10-4358	10-4368	5"	10-4380	10-4392
10-4394	10-4406	6"	10-4418	10-4430
10-4396	10-4408	6"	10-4420	10-4432
10-4398	10-4410	6"	10-4422	10-4434
10-4400	10-4412	6"	10-4424	10-4436
10-4402	10-4414	6"	10-4426	10-4438
10-4404	10-4416	6"	10-4428	10-4440


ACCESSORIES

PRODUCT NUMBER	MODEL	FEATURES
10-4478	Mid Size Wire Basket	16" x 10" x 1"
10-4480	Mid Size Wire Basket	16" x 10" x 2"
10-4482	Mid Size Wire Basket	16" x 10" x 3"
10-4484	Mid Size Wire Basket	16" x 10" x 4"
10-4975	Mid Size Silicon Mat	15" x 9"


GENERAL ACCESSORIES

Complete any SklarLite™ Container System by selecting disposable filters, indicators and labels to maximize security in preparation and maintenance of surgical sets and sterile materials.

PRODUCT NUMBER	MODEL	FEATURES
10-4900	Process Indicator Cards	Steam, 250/Bx
10-4910	Disposable Filters	100/Pk
10-4925	Security Seals	100/Pk
10-4930	ID Labels*	Silver
10-4930ENG	ID Labels*	Silver Custom Engraved, Max. 16 Characters

*Choose from six vibrant colors including Silver, Yellow, Red, Blue, Green or Black. Contact your Sklar Account Manager for information on the complete product line.



SKLARLITE™ FULL-SIZE CONTAINERS

23" L X 11" W

STANDARD CONTAINERS			SAFE CONTAINERS	
Non-Perforated Bottom	Perforated Bottom	Size	Non-Perforated Bottom	Perforated Bottom
10-4000	10-4012	4"	10-4024	10-4036
10-4002	10-4014	4"	10-4026	10-4038
10-4004	10-4016	4"	10-4028	10-4040
10-4006	10-4018	4"	10-4030	10-4042
10-4008	10-4020	4"	10-4032	10-4044
10-4010	10-4022	4"	10-4034	10-4046
10-4048	10-4060	5"	10-4072	10-4084
10-4050	10-4062	5"	10-4074	10-4086
10-4052	10-4064	5"	10-4076	10-4088
10-4054	10-4066	5"	10-4078	10-4090
10-4056	10-4068	5"	10-4080	10-4092
10-4058	10-4070	5"	10-4082	10-4094
10-4096	10-4108	6"	10-4120	10-4132
10-4098	10-4110	6"	10-4122	10-4134
10-4100	10-4112	6"	10-4124	10-4136
10-4102	10-4114	6"	10-4126	10-4138
10-4104	10-4116	6"	10-4128	10-4140
10-4106	10-4118	6"	10-4130	10-4142
10-4144	10-4156	8"	10-4168	10-4180
10-4146	10-4158	8"	10-4170	10-4182
10-4148	10-4160	8"	10-4172	10-4184
10-4150	10-4162	8"	10-4174	10-4186
10-4152	10-4164	8"	10-4176	10-4188
10-4154	10-4166	8"	10-4178	10-4190
10-4192	10-4204	10"	10-4216	10-4228
10-4194	10-4206	10"	10-4218	10-4230
10-4196	10-4208	10"	10-4220	10-4232
10-4198	10-4210	10"	10-4222	10-4234
10-4200	10-4212	10"	10-4224	10-4236
10-4202	10-4214	10"	10-4226	10-4238



ACCESSORIES

PRODUCT NUMBER	MODEL	FEATURES
10-4284	Full Size Wire Basket	21¼" x 10" x 1"
10-4286	Full Size Wire Basket	21¼" x 10" x 2"
10-4288	Full Size Wire Basket	21¼" x 10" x 3"
10-4290	Full Size Wire Basket	21¼" x 10" x 4"
10-4985	Full Size Silicon Mat	20½" x 9"

CANNULA CLEANING BRUSHES

PRODUCT NUMBER	FEATURES	QUANTITIES
10-1350	12" 2mm Diameter	3/pk
10-1660	12" 3mm Diameter	3/pk
10-1661	12" 3mm Diameter	12/pk
10-1663	12" 3mm Diameter	24/pk
10-1352	12" 4mm Diameter	3/pk
10-1354	12" 5mm Diameter	3/pk
10-1356	12" 6mm Diameter	3/pk
10-1664	15" 9mm Diameter	3/pk
10-1665	15" 9mm Diameter	12/pk
10-1666	15" 9mm Diameter	24/pk
10-1668	15" 25mm Diameter	3/pk
10-1669	15" 25mm Diameter	12/pk
10-1672	15" 25mm Diameter	24/pk
10-1362	18" 5mm Diameter	3/pk
10-1364	18" 10mm Diameter	3/pk
10-1366	18" 15mm Diameter	3/pk
78-2824	19" 8mm Diameter	3/pk
10-1370	24" 2.5mm Diameter	3/pk
10-1372	24" 5mm Diameter	3/pk
10-1374	24" 7mm Diameter	3/pk
10-1376	24" 10mm Diameter	3/pk



**INSTRUMENT CLEANING BRUSHES
(STAINLESS BRISTLES)**

PRODUCT NUMBER	FEATURES	QUANTITIES
10-1650	Stainless Steel Bristles	3/pk
10-1652	Stainless Steel Bristles	12/pk
10-1658	Stainless Steel Bristles	24/pk



**INSTRUMENT CLEANING BRUSHES
(NYLON BRISTLES)**

PRODUCT NUMBER	FEATURES	QUANTITIES
10-1651	Nylon Bristles	3/pk
10-1657	Nylon Bristles	12/pk
10-1659	Nylon Bristles	24/pk
10-1465	Nylon Bristles, Disposable	12/pk
10-1468	Nylon Bristles, Disposable	24/pk



**INSTRUMENT CLEANING BRUSHES
(DOUBLE ENDED)**

PRODUCT NUMBER	FEATURES	QUANTITIES
10-1444	Nylon Bristles	3/pk



10-1428 – ENDOSCOPIC BRUSH SET

PRODUCT NUMBER	DESCRIPTION	QUANTITIES
10-1350	Cannula Cleaning Brush 12" 2mm, 3/pk	1 pk
10-1660	Cannula Cleaning Brush 12" 3mm, 3/pk	1 pk
10-1354	Cannula Cleaning Brush 12" 5mm, 3/pk	1 pk
10-1362	Cannula Cleaning Brush 18" 5mm, 3/pk	1 pk
10-1364	Cannula Cleaning Brush 18" 10mm, 3/pk	1 pk
10-1370	Cannula Cleaning Brush 24" 2.5mm, 3/pk	1 pk
10-1372	Cannula Cleaning Brush 24" 5mm, 3/pk	1 pk
10-1376	Cannula Cleaning Brush 24" 10mm, 3/pk	1 pk
10-1398	Double End Valve Brush Nylon, 3/pk	2 pks
10-1444	Double End Brush Nylon, 3/pk	2 pks



SKLAR-CIDE™

HIGH LEVEL DISINFECTING/STERILIZING SOLUTION

Sporicidal • Bactericidal • Virucidal • Fungicidal

SKLAR-CIDE provides a safe and efficient means for the high-level disinfection or sterilization of medical, ophthalmic, dental and surgical devices including endoscopes and lensed instruments. It may be used as a standard soak solution or in automatic reprocessing units set at room temperature. Unlike other solutions, there is no need to increase the solution temperature to provide high-level disinfection in 25 minutes.

SKLAR-CIDE is a surfactant-free, activated glutaraldehyde liquid chemical sterilant and high-level disinfectant. Its active ingredient is glutaraldehyde in a nominal concentration of 3%. When the powdered activator is added, the solution turns green in color to indicate its readiness for use. In addition to the active ingredient, this ready-to-use solution contains buffers, wetting agents, and a rust and corrosion inhibitor and it has a pleasant fragrance.

SKLAR-CIDE produces high-level disinfection at room temperature in as little time as 25 minutes, and complete sterilization in 10 hours at room temperature.

- Reusable for 28 days instead of 14, safer and saves time and money.
- Surfactant free - Will not damage scope surfaces or cloud lenses.
- Can be used safely on all devices, including instruments and endoscopes. No need to stock several products.



PRODUCT NUMBER	FEATURES	QUANTITIES
10-1623	28 HLD 1 Gallon Bottles	4/case
10-1624	28 HLD 1 Gallon Bottle	1 bottle

GLUTARALDEHYDE STRIPS

PRODUCT NUMBER	FEATURES	QUANTITIES
96-4105	2.1% MEC 60 Strips / Bottle	2 bottles/ case



SKLAR DISINFECTANT

READY-TO-USE SPRAY SURFACE DISINFECTANT

- Ready to use, non-toxic formula.
- Specially formulated for the control of cross contamination.
- Effective against H1N1, HIV-1, HPV, VRE, Norovirus, Hep B & C as well as TB, MRSA, and A. Baumannii.
- Safe for all areas of the facility.
- Satisfies OSHA's Bloodborne Pathogen Standards.
- Non-irritating formula.
- EPA registered.

PRODUCT NUMBER	FEATURES	QUANTITIES
10-1643	24 oz. Spray Bottles	6/case
10-1642	24 oz. Spray Bottle	1 bottle
10-1653	1 Gallon Bottles	4/case
10-1654	1 Gallon Bottle	1 bottle
10-1644	2.5 Gallon Bottles	2/case



SKLAR DISINFECTANT SURFACE WIPES

READY-TO-USE, PRE-MOISTENED TOWELETTES

**SAME EFFECTIVE FORMULA AS
SKLAR DISINFECTANT!**

- Ready to use, non-toxic formula.
- Specially formulated for the control of cross contamination.
- Effective against H1N1, HIV-1, HPV, VRE, Norovirus, Hep B & C as well as TB, MRSA, and A. Baumannii.
- Safe for all areas of the facility.
- Satisfies OSHA's Bloodborne Pathogen Standards.
- Destroys odor-causing bacteria.
- 7" x 10" non-woven, tear resistant towelette.
- Clean scent.
- EPA registered.



KILLS MRSA AND H1N1 VIRUS

PRODUCT NUMBER	FEATURES	QUANTITIES
10-1616	100 Towelettes / Canister	6/case
10-1617	100 Towelettes / Canister	12/case

SKLAR SPRAY-ZYME™

SURGICAL INSTRUMENT PRE-CLEANER

Multi-Enzyme Based Instrument Presoak & Cleaner Foam Spray Cleaner

The newest product from Sklar for expert care of your surgical instruments. This multi-enzymatic spray creates a heavy foam which lasts for hours, keeping soiled instruments moist to avoid encrustation while also pre-cleaning prior to processing.

- Prevents Organic Encrustation - enzyme foam stops blood and bodily fluids from drying onto surgical instruments.
- No More Spills - messy spills between the operating room and central cleaning /processing department are eliminated.
- Neutral pH - safe and effective for stainless steel, tungsten carbide, aluminum and brass instruments.
- Bacteriostatic - microbiological control agents inhibit bacterial growth in the ready-to-use foam solution.
- Long Lasting Foam - saves time, foam lasts many hours after application.
- Multi-Enzymatic - contains a synergistic blend of lipase and protease enzymes.

PRODUCT NUMBER	FEATURES	QUANTITIES
10-2722	22 oz. Refillable Spray Bottles	12/case
10-2724	1 Gallon Bottles	4/case

SKLAR-ZYME™ POUCHES

CONVENIENT VALUE

Each pouch of SKLAR-ZYME produces one gallon of enzymatic detergent. This convenient, easy to use product guarantees very little waste and is ideal for the office or smaller surgical setting.

PRODUCT NUMBER	FEATURES	QUANTITIES
10-2774	½ oz. Pouches	32/box



SKLAR SOAK™

INSTRUMENT PRE-SOAK & DECONTAMINANT RUST RESISTANT!

Immersing soiled instruments in SKLAR SOAK effectively pre-cleans and decontaminates, without fear of rusting!

SKLAR SOAK is an effective pre-cleaner that breaks down bio-burden, lifting blood and other organic debris from the surface, while penetrating box locks and crevices. The rust inhibitor protects the delicate chromium oxide layer of the instrument, thus defending against rust and corrosion. Using SKLAR SOAK as a pre-cleaner is the first line of defense against rust, prolonging the life of the instrument and therefore protecting your investment!

Soaking for only two minutes thoroughly decontaminates instruments, reducing the risk of exposure to pathogens such as H1N1, HIV-1, HPV, VRE, Norovirus, Hep B & C as well as TB, MRSA and A. Baumannii. However, because SKLAR SOAK contains a rust inhibitor, instruments may soak for several hours without fear of rusting! Employees are protected from cross contamination, and at the same time can enjoy the freedom to reprocess instruments when it's convenient!

SKLAR SOAK is ideal for use on all stainless steel instruments:

- PH balanced, designed specifically for stainless steel instruments.
- Intermediate level disinfectant with broad spectrum of efficacy.
- Rust resistant.
- Non-caustic.
- Ready to use.
- No special handling or disposal.
- Fragrant, no offensive odors.
- Non-glutaraldehyde.



PRODUCT NUMBER	FEATURES	QUANTITIES
10-1603	1 Gallon Bottles	4/case
10-1609	1 Gallon Bottle	1 bottle

Sklar is pleased to offer a full line of Care & Cleaning agents to disinfect, decontaminate, sterilize, lubricate and protect your stainless steel surgical instruments and equipment. For expert care of your surgical instruments at every stage of the cleaning process, trust Sklar to deliver the right cleaning agent for you.



3 EASY STEPS FOR THE CARE AND CLEANING OF YOUR INSTRUMENTS

1 DECONTAMINATION

Sklar Disinfectants are a pre-cleaner and disinfectant for the safe decontamination of soiled surgical instruments, medical and surgical equipment and inanimate surfaces including plastics, stainless steel, chrome, glass, Plexiglas and baked on painted surfaces.

Disinfectants are effective against a broad range of pathogens in as little as 2 minutes.



2 CLEANING

Sklar carries a complete line of effective detergents and enzymatics for cleaning, not only for your surgical instruments, but also metals, glassware, plastic and rubber. Available in your choice of powder, liquid concentrate or spray foams.



3 LUBRICATION

To guarantee extended life and performance of your surgical instruments, choose one of Sklar's lubricants or sprays as the final step in your instrument processing regimen. Specifically designed to protect against corrosion, rust and staining, Sklar offers several different types of lubricants to protect your investments.





889 South Matlack Street
West Chester, PA 19382

ph: 800.221.2166
610.430.3200
fx: 610.696.9007
e-mail: surgi@sklarcorp.com
e-mail: international@sklarcorp.com

www.sklarcorp.com



DIGITAL FILE ONLY